

MARATHON LOGGER PULSE™

...an Enterprise Software Solution for Managing The Cold Chain



Our enterprise software is the ideal solution for managing and analyzing data from multiple data loggers from multiple locations.

NO LIMITS – Any location. Any time. Any users.

Features:

Location Independence: Data loggers located anywhere can be centrally managed and monitored.

Flexible Hosting Options: From single Windows PC to multiple servers in multiple locations; local in-house servers to cloud-based services.

Web-base Software: Means hassle-free deployment for access by an unlimited number of users with a standard web browser.

Dashboard Control: An intuitive and interactive dashboard makes it easy to view status, manage, and program multiple loggers with ease.

Alarm Notifications: Highly customizable alarm notifications ensure that the right people are alerted when loggers detect a problem. Alarm notifications can be sent as

21 CFR COMPLIANT

- Controls for closed and open systems
- Electronic audit trails for regulatory compliance.
- Electronic audit trails for regulatory compliance.

Name	Last Mean	Last Mean Date	Configuration	Last Status Date	Notes
edNET 102 F 804804	20.2	4/18/2012 9:21:32 AM	30 sec manual alarm 10-d	4/18/2012 9:22:30 AM	(179)
edNET 103	18.7	4/18/2012 9:21:14 AM	1 min Internal	4/18/2012 9:22:40 AM	(8)
edNET 104 F 804804	19.8	4/18/2012 9:20:04 AM	1 Min & Alarm	4/18/2012 9:22:20 AM	(150)
edNET 106 F 804804	19.1	4/18/2012 9:19:26 AM	1 Min & Alarm	4/18/2012 9:22:29 AM	(170)
edNET 107 F 804804	21.2	4/18/2012 9:21:18 AM	1 min Internal	4/18/2012 9:22:41 AM	(16)
edNET 109 F 804804	16.5	4/18/2012 9:21:10 AM	1 Min & Alarm	4/18/2012 9:22:30 AM	(17)
edNet 110	17.9	4/18/2012 9:20:44 AM	1 Min & Alarm	4/18/2012 9:22:40 AM	(8)
edNET 111	19.1	4/18/2012 9:14:58 AM	1 min Internal	4/18/2012 9:22:28 AM	(7)
edNET 112	20.8	4/18/2012 9:18:58 AM	1 min Internal	4/18/2012 9:22:28 AM	(7)
edNET 113	18.9	4/18/2012 9:18:58 AM	1 min Internal	4/18/2012 9:22:28 AM	(7)
edNET 114	21.6	4/18/2012 9:20:12 AM	1 min Internal	4/18/2012 9:22:29 AM	(8)
edNET 115	21.5	4/18/2012 9:18:34 AM	1 min Internal	4/18/2012 9:22:28 AM	(7)
edNET 116	21.1	4/18/2012 8:58:14 AM	1 min Internal	4/18/2012 9:22:28 AM	(7)
edNET 117	20.6	4/18/2012 9:15:14 AM	1 min Internal	4/18/2012 9:20:19 AM	(1)
edNET 119	19.1	4/18/2012 9:20:34 AM	1 min Internal	4/18/2012 9:22:28 AM	(7)
edNET 121	21.2	4/18/2012 9:20:22 AM	1 Min & Alarm	4/18/2012 9:22:39 AM	(7)
edNET 122	20.8	4/18/2012 9:20:56 AM	1 min Internal	4/18/2012 9:22:40 AM	(2)
RF4866	19.7	4/18/2012 9:22:16 AM	1 min Internal	4/18/2012 9:22:19 AM	(6)
RF4867	19.2	4/18/2012 9:22:16 AM	1 min Internal	4/18/2012 9:22:16 AM	(6)
RF4871	20.3	4/18/2012 11:02:16 AM	1 min Internal	4/18/2012 9:20:42 AM	(8)
RF4874	22.1	4/18/2012 9:16:58 AM	1 min Internal	4/18/2012 9:22:28 AM	(8)
RF4881	-19.4	4/18/2012 9:19:02 AM	3 Min	4/18/2012 9:22:34 AM	(13)

Dashboard Control

a phone call to land lines and mobile phones, as email to regular email accounts or to mobile phones as text messages.

Real-time Analytics: Multiple users can concurrently monitor data streamed in real-time from any number of loggers.

Real-time Data: Logger data can be displayed individually or aggregated.

Bulk Administration: Programming, starting and stopping hundreds of loggers can be done effortlessly from a single web page.

Concurrent & Secure Access: Enable multiple users to simultaneously perform analytics and administration functions while constraining actions based on permissions.


Audit Trail: Individual actions and exceptions are logged to the database for accountability and traceability. User notes are incorporated for explanation and follow up.

Historical Analytics: View aggregated logger data for any period in the past in graphical or tabular formats. Export to Microsoft Excel for further analysis and reporting.

Attribute Filters: A single click is all it takes to filter for specific status and conditions, e.g. Alarms.

Multi-Sensor Inputs: Temperature, Humidity, Vibration Monitoring, Electrical Power Consumption, CO2, N2, O2 (Pending).

Support for all popular browsers: Flexible, intuitive interface ensures that LoggerPulse will work the same way regardless of your business browser standard.

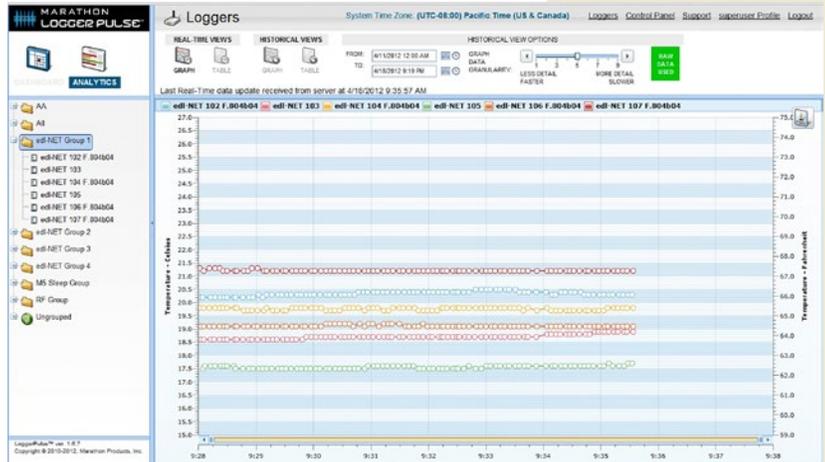

Real-Time Analytics



Historical Analytics

Please call 1-510-562-6450
for more information, a demonstration, or to place an order.

www.marathonproducts.com

627 McCormick St., San Leandro, CA 94577 ■ sales@marathonproducts.com
Specifications subject to change without prior notice.